



SINCE 1966
MADE IN ITALY

Wi.Neb

by **FLAEM**®



**Powerful and compact
aerosol therapy unit**

Speed, Silence, Comfort

**New RF7 Dual Speed Plus
nebulizer**

Use with Mouthpiece




Maximum therapeutic effect is ensured by the valve system allowing optimisation of medication delivery and particle size to reach the lower respiratory tract.

Minimum dispersion of medication is ensured by the synchronised opening and closing of the valves. Exhaled air is directed towards the dedicated exit. The valve limits dispersion and the aerosol mixture is enriched for the next inhalation.

Air filtration

Air filtration is ensured by the special filter. Easy to reach, it is located under the little lid on the side that together with the casing acts a **silencer**.

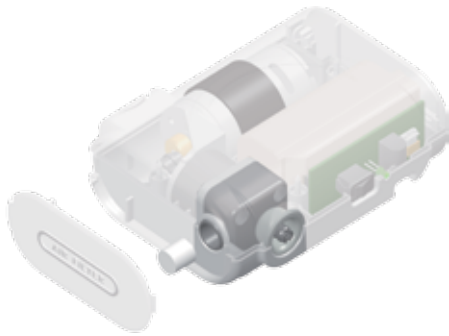

SINCE 1966
MADE IN ITALY
Wi.Neb
by **FLAEM**®



Why choose Wi.Neb



RF7 Dual Speed Plus nebulizer



high performance



extremely silent



adult/child SoftTouch masks



environmentally friendly
low energy consumption



usable in automobiles




convenient carrying bag

Wi.Neb with power supply unit or cigarette lighter cable

Compact size and powerful device. It can be powered by the Lithium Power battery or by a 100/240V - 50-60Hz multi voltage power supply unit or by a cigarette lighter power cord.

Between 2 and 6 μm
for laryngotracheal
deposition

Between 0.5 and 3 μm
for alveolar deposition

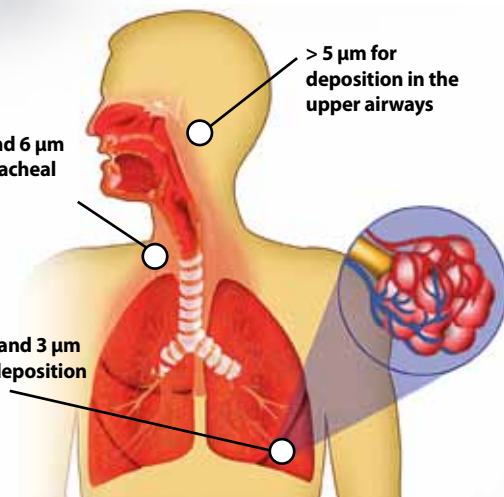

MMAD (1): Mass Median Aerodynamic Diameter
FPF (1): Fine Particle Fraction (breathable fraction)
ml/min (2): delivery speed

(1) In vitro testing certified by
TÜV Rheinland LGA Products GmbH - Germany
in compliance with the **European Standard EN 13544-1**
for nebulizing systems

(2) Data shown is in accordance with **Flaem Nuova** internal
procedure **I29-P07.5**. Values shown for the aerosol delivery
may vary depending on the patients respiratory capacity

In compliance with:
European Standard **EN 10993-1**
"Biological evaluation of medical device"
and **93/42/EEC Medical Device Directive**



FLAEM NUOVA S.p.A.
Company with certified quality system



Export Department



medexpo international s.r.l.

Fax: (02) 70991111